


EZEEKLEEN 2.5

Distillation Alternative

Optional Ultraviolet Cartridge

**Turn your
Ezee Kleen 2.5
into an Ultra
Pure Water
System!**

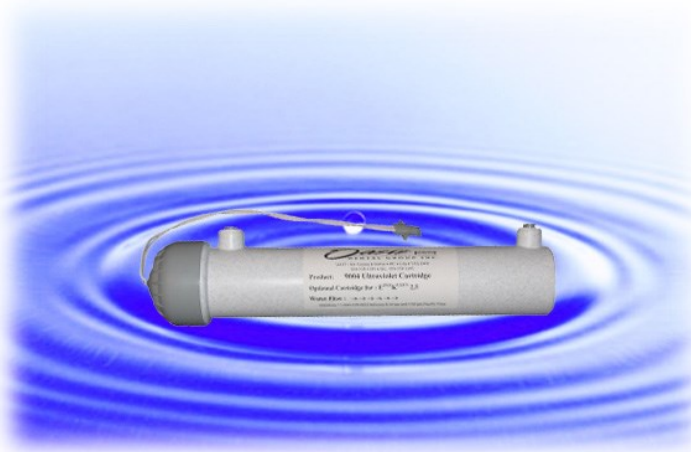

**Produce water
Ideal for your Self
Contained Dental
Unit Water
Systems**

Why Use Ezee Kleen 2.5 Water with our Optional Ultraviolet Cartridge in your Self Contained Dental Unit Water System?

Destroy Microbial Contamination Before it Becomes a Problem!

The Ezee Kleen 2.5 Ultraviolet (UV) Cartridge produces 32,000 w-sec/cm² of UV radiation which destroys bacteria, viruses and other micro-organisms by interfering with the DNA and RNA in the organisms' reproductive cycles. Effective on bacteria, algae, viruses & yeasts.

Reduce the change frequency of micro-filtration cartridges on your dental lines to save time and money.

No organic matter (algae, etc.) **and inorganic elements** (dissolved solids - sodium, calcium, etc.) means no food and nutrients for bacteria survival and biofilm development. Slower growth if it does occur.

Reduce chemical treatment frequency on your dental lines to save time, money and deterioration of lines, o-rings, gaskets and equipment from exposure to harsh chemicals.

No Microbes means less chance of biofilm growth in the first place

Additional Benefits of Quality Ezee Kleen 2.5 Water!

Autoclave Water - onsite, on demand - 1 gallon (4 litres) in 2 minutes of consistently monitored water to ensure quality complies with all autoclave manufacturer's TDS level guidelines

Eliminate the bulky storage of bottles

Ensure you never run out of water for your autoclaves - Eliminate last minute store runs for water.

Self monitor for consistent water quality (Unlike bottled water)

Maximize the efficiency & life of dental unit equipment and instruments. Quality water ensures your equipment does not become clogged with mineral scale that can cause malfunction, increased energy usage and shortened life.

Save Money!!!

Fast Simple Installation

Use the 9004 Ultraviolet Cartridge Install Kit to install the Ultraviolet Cartridge for the first time. Takes only minutes to install. No special tools required.

9004 Ultraviolet Cartridge Install Kit (Includes one lamp)

9004-01 Ultraviolet Replacement Lamp (Change Annually)



15531 - 24th Avenue, Unit 24, Surrey, BC V4A 2J3 • Ph: (800) 338-6693 • (604) 538-1285 • Fax: (604) 538-1295

ORGANISMS DESTROYED BY Ezee Kleen 2.5 ULTRAVIOLET RADIATION CARTRIDGE GERMICIDAL LEVEL (w-sec/cm²)

The **EZEKLEEN 2.5 Ultraviolet Cartridge** effectively uses Ultraviolet (UV) radiation to destroy bacteria, viruses and other microorganisms by interfering with the DNA and RNA in the organisms' reproductive cycles. The contact disinfectant property of UV radiation energy is measured in microwatt-seconds per square centimetre (w-sec/cm²).

BACTERIA	BACTERIA, continued
Agrobacterium8,500	Shigella dysenteriae (Dysentery)4,200
Bacillus anthracis.....8,700	Shigella flexneri (Dysentery).....3,400
Bacillus megaterium (vegetative)2,500	Shigella sonnei.....7,000
Bacillus subtilis (vegetative).....11,000	Staphylococcus epidermidis.....5,800
Clostridium tetani22,000	Staphylococcus aureus7,000
Corynebacterium diphtheriae.....6,500	Streptococcus faecalis.....10,000
Escherichia coli.....7,000	Streptococcus hemolyticus5,500
Legionella bozemanii.....3,500	Streptococcus lactis.....8,800
Legionella dumoffii5,500	Viridans streptococci3,800
Legionella gormanii.....4,900	Vibrio cholerae6,500
Legionella micdadei.....3,100	
Legionella longbeachae2,900	ALGAE
Legionella pneumophila3,800	Chlorella vulgaris (algae).....22,000
Leptospira interrogans (Infectious Jaundice)6,000	
Mycobacterium tuberculosis.....10,000	VIRUSES
Neisseria catarrhalis.....8,500	Bacteriophage (E.coli)6,600
Proteus vulgaris6,600	Hepatitis virus8,000
Pseudomonas aeruginosa laboratory strain3,900	Influenza virus6,600
Pseudomonas aeruginosa environmental strain10,500	Poliovirus (Poliomyelitis)21,000
Rhodospirillum rubrum6,200	Rotavirus.....24,000
Salmonella enteritidis7,600	
Salmonella paratyphi (Enteric Fever)6,100	YEAST
Salmonella typhimurium.....15,200	Baxter's yeast.....8,800
Salmonella typhosa (Typhoid Fever).....6,000	Brewer's yeast6,600
Sarcinia lutea26,400	Common cake yeast.....13,200
Serratia marcescens6,200	Saccharomyces var. ellipsoideus13,200
	Saccaromyces sp.....17,600

EZEKLEEN 2.5 UV Cartridge

Includes all componentry necessary for initial installation, includes 9004-01

Product ID: 9004

Dimensions: 11.5" x 2.25"/ 29cm x 6 cm

Weight: 1 lb/0.45 kg

Replacement Lamp: 9004-01

Recommended change frequency: 9000 hrs or annually

Power Supply: 55 VAC; 60 Hz; 16 amp. (Optional European format available)

Approved Componentry: UL/CSA approved, made of FDA approved materials and requires no MSDS reporting


These data are believed to be reliable but do not imply any warranty or performance guarantee. Actual performance will vary.
Microbial testing may be appropriate to determine suitability of water for intended use.

FZEEKLEEN 2.5

Distillation Alternative

Model: 9000

OPTIONAL 9004 UV CARTRIDGE INSTALLATION DIAGRAM

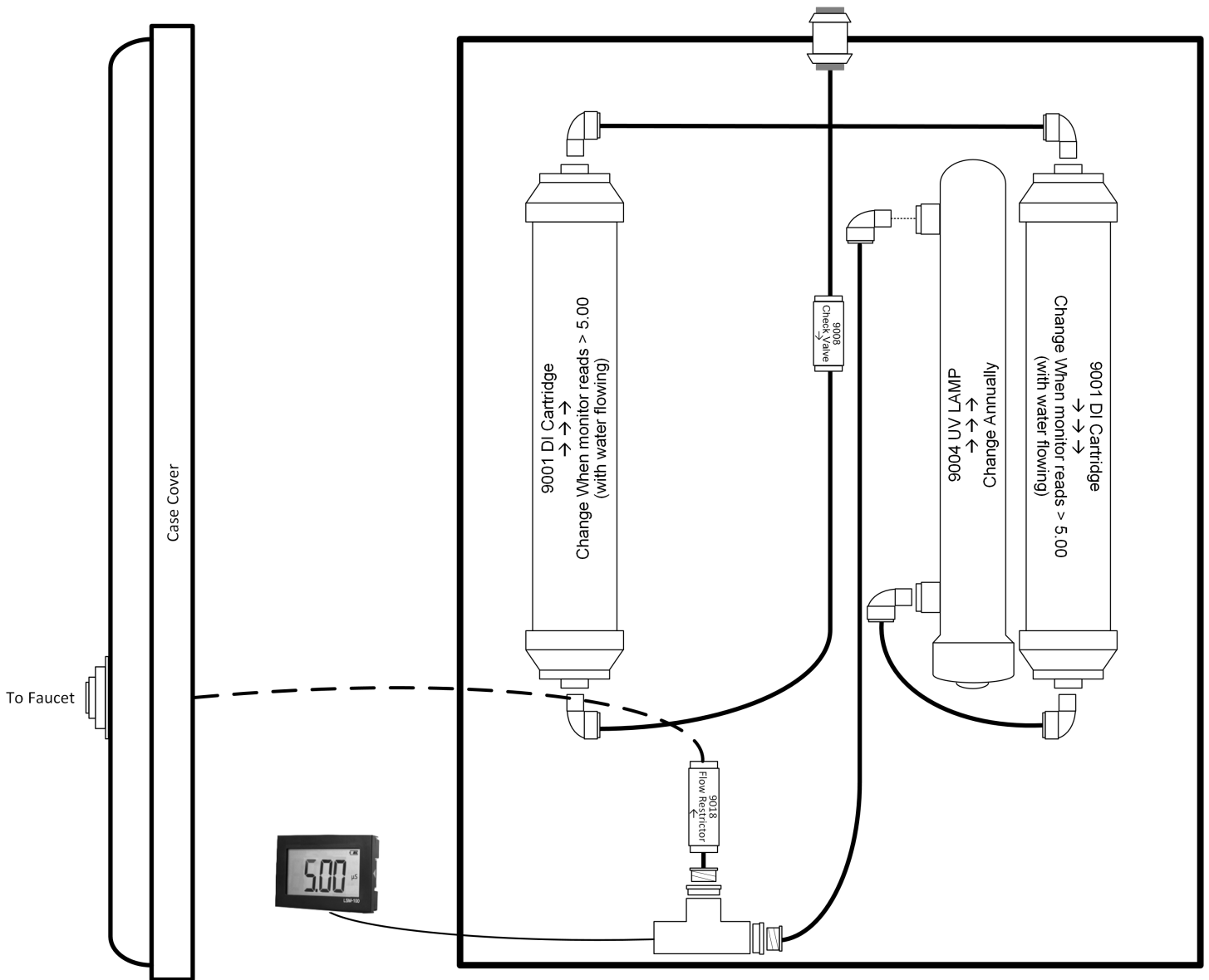

Questions? Contact customer support Mon – Fri, PST, 8:00am – 5:00pm.



15531 - 24th Avenue, Unit 24, Surrey, BC V4A 2J3 • Ph: (800) 338-6693 • (604) 538-1285 • Fax: (604) 538-1295

© 2009-2024 Oasis Dental Group Inc.
www.oasisdentalgroup.ca